

## SANITIZED PRESSURE RELIEF CUSHIONS



### Description

Sanitized suapel anti-decubitus cushions filled with conjugated silicone fibre. With straps to fix the cushion to the wheelchair or other seats. Available in round with hole, square, square with hole, square horseshoe and round horseshoe.

### Effects and advantages

Thanks to its composition and shape, the pressure area is relieved, preventing the appearance of pressure ulcers (UPP).

Recommended for people who remain seated for long periods of time due to reduced mobility, sedentary life, etc. Also, depending on the pathology (postpartum, haemorrhoids, fistulas, prostate operations, etc.), one form or another is recommended to relieve pressure on the affected area.

Sanitised treatment of the tissue:

Anti-mite, anti-bacterial, anti-fungal, anti-bugs and anti-odour.

Maximum support weight: 80kg.

Risk level: Low-medium.

### Installation

The cushion fits any type of seat, but should be checked first. Please note that for horseshoe-shaped cushions, the open part is positioned in the backrest area (see illustration).

### Composition

- Outer cover:

Suapel fabric (92% polyester and 8% acrylic) with sanitized treatment. Finished with acrylic resins and antibacterial sanitized.

Inner part: terlenka fabric (65% polyester and 35% cotton).

- Filling: Conjugated siliconised hollow fibre.

It has 100% polyester straps to attach the cushion to the wheelchair or other seats and facilitate its position.

### Dimensions

- Round with hole: 44x11cm.

- Other cushions: 44x44x11cm.

### Maintenance

Periodically brush or shake the fibre face to keep it in optimal condition and slow down deterioration. Hand wash with cold water and non-abrasive products.

In case of machine washing, use a programme for delicate garments and do not spin-dry.

It is recommended to use a special cover or net to avoid excessive deformation during washing.

Dry flat.

Do not iron.

Do not bleach.

Do not dry clean.


### Warranty

The warranty period for all cushions except Liquid Flotation and Silicone cushions is two years from the date of purchase of the cushion. The Liquid Flotation and Silicone cushions are guaranteed for 3 months. This warranty does not cover normal wear and tear of the product and does not replace legal warranties. Please contact the company's quality department in the event of an incident in order to resolve it.


# Anti-decubitus

## Sanitized pressure relief cushions


### Models

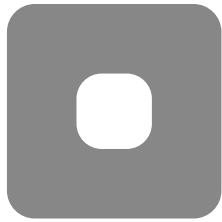

Square  
w/hole

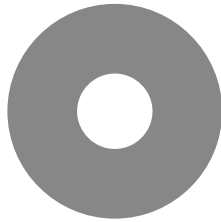


Round  
w/hole

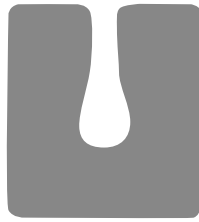

Square  
horseshoe



Round  
horseshoe

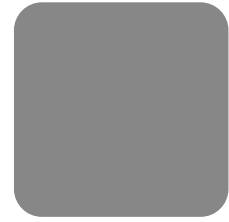

Square

### Colours



Navy  
blue



Grey



White

Correct way to  
place a horseshoe  
cushion.



MODEL	DIMENSIONS (CM)	COLOUR	VERSION	REF.
Sanitized pressure relief cushions	44x11	White	Round w/hole	101100
	44x44x11	White	Square w/hole	101101
			Round horseshoe	101102
			Square	101103
			Square horseshoe	101104
			Grey	Square
	44x11	Navy blue	Square	101108
		Grey	Round w/hole	101200
	44x44x11	Navy blue	Round w/hole	101201
		Grey	Round horseshoe	101202
		Navy blue	Round horseshoe	101203
		Grey	Square horseshoe	101204
		Navy blue	Square horseshoe	101205
		Grey	Round w/hole	101210
Navy blue		Round w/hole	101212	



antiol or



anti bacterias



anti hon gos



anti ácar os



anti CHINCHES



SUAVE



TRANSPIRABLE



LAVABLE